

بسم الله الرحمن الرحيم
جامعة أم درمان الإسلامية
كلية الصيدلة
قسم علم الأدوية
الفرقة الأولى
معمل أحياء الخلية والوراثة

Lab NO (1)

Lab Safety in Biological Sciences Laboratories

Laboratory Safety

Safety in the Laboratory Should always be in your mind. Throughout this manual Safety recommendations are given, below are some general consideration that anyone in a laboratory should know.

▪ General laboratory safety precaution.


1. Follow all instructions carefully given by your teacher.
2. Determine the location of Fire Extinguishers, Chemical safety showers and Eye washers, Chemical Spill Kits, and alternative exit routes for lab evacuation.
3. Remember that eating, or drinking in the lab room is totally prohibited.
4. Wear lab aprons when working with chemicals, hot material, or preserved specimens.
5. Wear safety goggles when using dangerous chemicals, hot liquids, or burners.
6. Any chemicals spilled on the hands or other parts of the skin should be washed off immediately with a plenty of running water.
7. If you have an open skin wound, be sure that it is covered with a water proof bandage.
8. Never work alone in the laboratory.
9. Keep your work area clean & dry.
10. Turn off all electrical equipment, water, and gas when it is not in use, especially at the end of the laboratory period.

▪ Special precautions for working with heat or fire:

1. Never leave a lighted Bunsen burner or hot object unattended. When an object is removed from the heat & left to cool, it should be placed where it is shielded from contact.



Microscope Parts

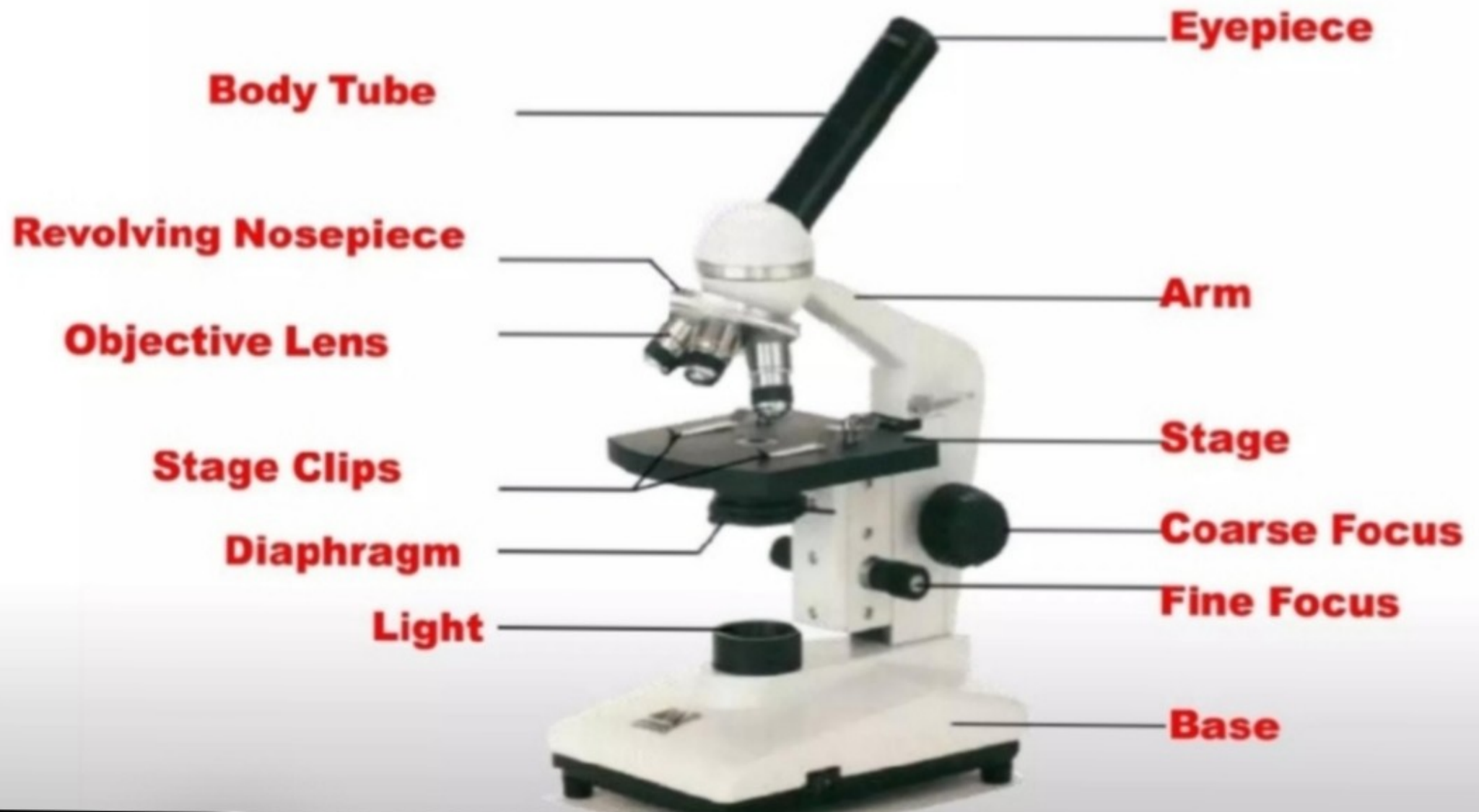

- **Eyepiece Lens:** the lens at the top that you look through.
- **Body Tube:** Connects the eyepiece to the objective lenses.
- **Objective Lenses:** Usually you will find 3 or 4 objective lenses on a microscope. They almost always consist of 4X, 10X, 40X and 100X powers.
- **Nosepiece:** This is the part that holds the objective lenses and can be rotated to change power.
- **Base:** The bottom of the microscope, used for support.
- **Arm:** Supports the tube and connects it to the base

Light source: found near the base of the microscope; makes the specimen easier to see.

Stage: The flat platform where you place your slides.

Stage clips: found on the stage; hold the slides in place.


Diaphragm: rotating disk under the stage; different sized holes vary the intensity and size of light that is projected upward into the slide.

Coarse adjustment knob: large, round knob on the side of the microscope used for focusing the specimen; it moves the stage up and down.

Fine adjustment knob: small, round knob on the side of the microscope used to fine-tune the focus of your specimen after using the coarse adjustment knob.




٢٠ نوفمبر، ٢٠٢٥
٧:١٣ م →

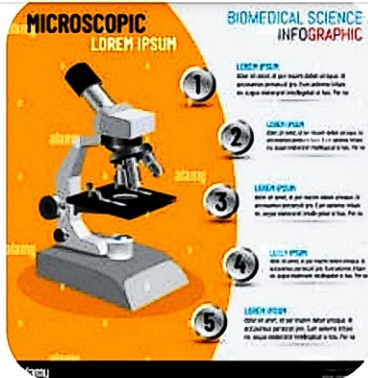

To use a
microscope,
start with the
lowest power
objective lens

and place the slide on the
stage, centering the specimen.
Then, while looking through
the eyepiece, use the coarse
focus knob to bring the image
into view, followed by the fine
focus knob for a sharp image.
Finally, center the desired area,
and rotate to higher power
objectives, readjusting with
the fine focus as needed.

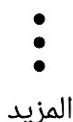

المزيد



حذف



تعديل



المفضلة



مشاركة

